

INDIBA®
revitalizing lives



REHABILITATION & SPORTS

INTIMATE HEALTH

COMPACT PRO COMPACT PLUS

—
Are You Ready to Play?


www.indiba.com

COMPACT Plus & Pro

Power, Precision, and Versatility

Advanced radiofrequency devices designed for premium rehabilitation, sports, intimate health, and dermatofunctional treatments.

Built for professionals who demand high performance, these devices offer exceptional power, precision, and versatility to deliver superior clinical results.

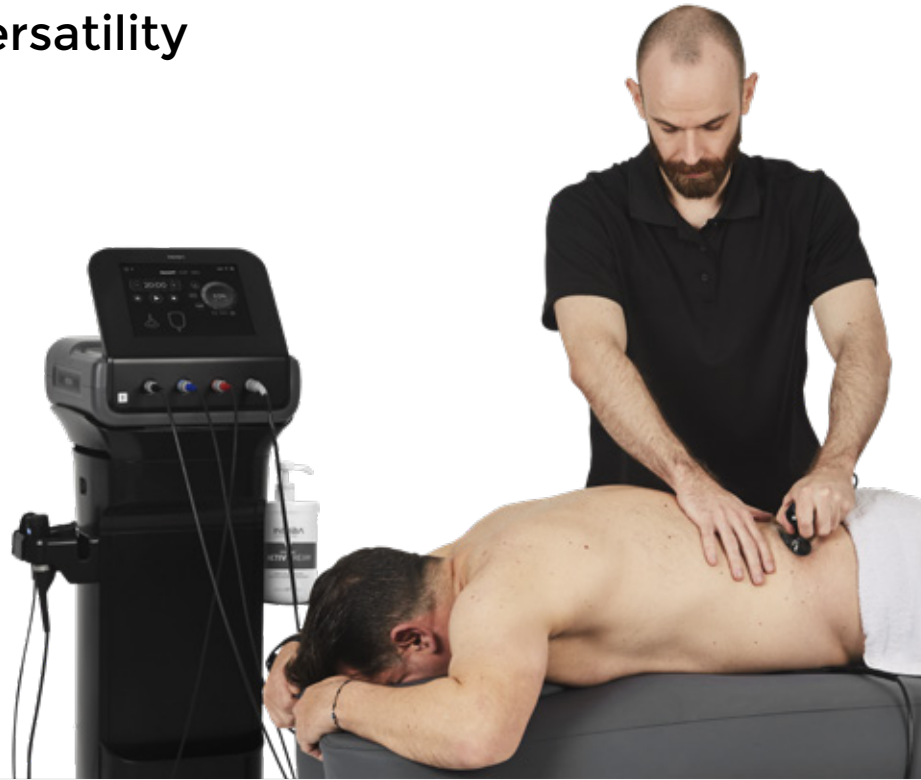

Exclusive DUAL RES Mode (COMPACT PRO only)

Treat multiple areas simultaneously using two RES electrodes with one return plate.

Real-Time Temperature Control

LED feedback and handle buttons provide instant, at-hand temperature monitoring and power adjustments.

Smart Handle for Seamless Treatments

Single applicator for CAP and RES, enabling quick electrode changes.



COMPACT PRO

The Science Behind INDIBA's Technology

Choosing INDIBA device means investing in the gold standard of patient care—trusted, reliable, and scientifically proven for superior results. Its stable 448 kHz frequency ensures consistent, clinically validated outcomes, backed by research for predictable, effective treatments.

Proven Effects of INDIBA's 448 kHz Technology:

Thermal Effects

- Increases metabolism
- Improves circulation and drainage
- Stimulates collagen & elastin synthesis
- Relieves muscle tension and pain

Subthermal Effects

- Promotes stem cell proliferation & differentiation
- Enhances tissue repair
- Reduces inflammation and fibrosis
- Increases metabolism
- Regulates cell migration



High Power for Optimized Results

Effectively treats chronic conditions and large areas with consistent outcomes.

Smart Electrodes for Precision Therapy

Built-in temperature sensors ensure controlled, accurate, and comfortable treatments.

Comprehensive Accessory Range

Enhance treatments with specialized accessories for rehab, sports, intimate health, and dermatofunctional therapy.

COMPACT PLUS



COMPACT PLUS



COMPACT PRO

Device

Power Output

Average Power: 120 W RES / 375 VA CAP
Peak power: 240 W RES

Average Power: 180 W RES / 425 VA CAP
Peak power: 360 W RES

Smart Technology

Smart handle & temperature sensors

Smart handle, temperature sensors, Dual RES mode

Weight

7.3 Kg

7.3 Kg

Touchscreen Size

10.4-inch capacitive touchscreen

10.4-inch capacitive touchscreen

Dual RES Mode

No

Yes

Temperature Sensors

Yes

Yes

Smart Handle

Yes

Yes

Software Features

Advanced software with real-time monitoring, treatment guides

Sophisticated software with real-time monitoring, treatment guides, temperature visual scale and DUAL RES mode.

Accessories Included

1 Smart Handle (RES/CAP) with cable
1 Return Plate
3 Smart CAP Electrodes (Ø 40, 55, 65 mm)
3 Smart RES Electrodes (Ø 35, 50, 65 mm)
Trolley

1 Smart Handle (RES/CAP) with cable
1 Return Plate
4 Smart CAP Electrodes (Ø 40, 55, 65, 80 mm)
3 Smart RES Electrodes (Ø 35, 50, 65 mm)
Trolley
DUAL RES:
→ 25 Adhesive RES return plates,
→ Smart cable, RES cable, Return cable

Add-on Accessories Available

Mobility Pack:
→ 2 Small CAP & RES Handles with cables
→ 2 CAP electrodes (Ø 40 & 60 mm)
→ 2 RES electrodes (Ø 35 & 50 mm)
Essence kit
Pluma Smart
Longitudinal CAP and RES electrodes
Cylinder return plate
New Fascia pack (Max & Flex) with magnetic cable
Hands-Free electrodes
Additional Smart electrodes:
→ 2 CAP Smart electrodes (Ø 19, 30 mm)
→ 1 Smart RES electrode (Ø 19 mm)

Mobility Pack:
→ 2 Small CAP & RES Handles with cables
→ 2 CAP electrodes (Ø 40 & 60 mm)
→ 2 RES electrodes (Ø 35 & 50 mm)
Essence kit
Pluma Smart
Longitudinal CAP and RES electrodes
Cylinder return plate
New Fascia pack (Max & Flex) with magnetic cable
Hands-Free electrodes
Additional Smart electrodes:
→ 2 CAP Smart electrodes (Ø 19, 30 mm)
→ 2 Smart RES electrodes (Ø 19, 90 mm)

Exclusive Features for Personalized Treatments



DUAL RES Mode

The DUAL RES Mode alternates two RES electrodes with one return plate for efficient, dual-area treatment. Designed for rehab, sports recovery, and injury treatment, it enhances therapy with Active RF, targets complementary muscles, and enables bilateral, time-saving treatments. For Intimate treatments, it supports core stability, posture, and pelvic recovery.




Smart Handle

- Dual-Function: Switches seamlessly between CAP & RES treatments.
- Real-Time Temperature Control: Ensures precision and safety.
- Power Control At-Hand: Handle-integrated power control allows intensity changes without interruption, keeping the focus on the patient.
- Ergonomic & Durable: Built for long-term comfort.



New Fascia Pack (Max & Flex) with Magnetic Cable

Our stainless-steel fascia electrodes provide exceptional versatility, durability, and precision. Designed to effectively treat myofascial restrictions, they enhance safety, performance, and ease of use.

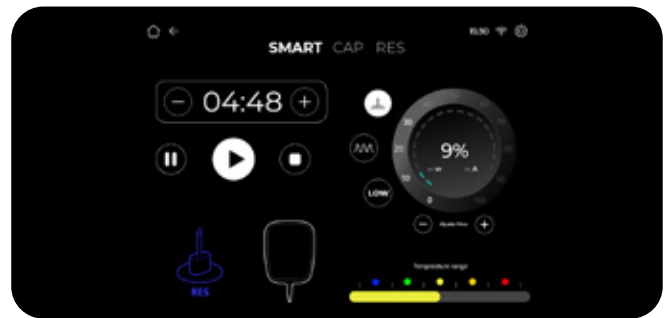

Precision & Control for Seamless Treatments

LED color feedback on the handle and a visual temperature scale on screen provide instant temperature monitoring for accurate and safe treatments.

Smart Electrodes

- Temperature Monitoring: Ensures accurate, safe treatments.
- Enhanced Comfort: Optimized for effective patient care.



Applications in Rehabilitation & Sports


INDIBA technology delivers targeted treatment for pain relief, injury recovery, and musculoskeletal conditions:



- **Chronic Pain:** Eases arthritis and joint discomfort.
- **Acute Injuries:** Speeds up recovery from sprains and strains.
- **Musculoskeletal Disorders:** Treats cervical, lower back pain, and more.
- **Post-Surgical Care:** Promotes tissue healing and function restoration.
- **Tendinopathy:** Reduces pain and inflammation.

Specialized Accessories

- **Hands-Free Electrodes:** Enables contact-free radiofrequency therapy for the user, while allowing active treatment and exercise for the patient.
- **Precision Electrode Pack:** Optimized for small, contoured, and hard-to-reach areas.
- **CAP & RES Additional Electrodes:** Customizable sizes and shapes ensuring precise and effective treatments.

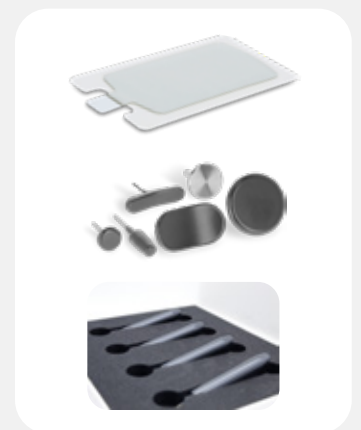

INDIBA Creams: Essential for Optimal Performance

Designed and tested to ensure excellent manoeuvrability of the electrodes and the effectiveness and safety of the treatment.

- Maximizes treatment results.
- Enables smooth application and handling.
- Ensures safe, reliable performance.
- Preserves device warranty for long-term protection.



Intimate Health and Dermatofunctional



Intimate Health

- Pelvic pain conditions
- Vaginal atrophy & rejuvenation*
- Postpartum recovery
- Post-surgical pain & edema
- Improved quality of life

*Not under Notified Body certification. These indications are not in the scope of medical device regulation (EU) 2017/745.



Dermatofunctional

- Skin injury healing
- Skin rejuvenation* (wrinkle reduction, skin firming and double chin reduction)
- Improvement of the appearance of orange-peel skin
- Lymphatic and venous drainage

Specialized Solutions for Intimate Health



Pluma Smart Intracavitary Handle

Advanced temperature-controlled device with single-use CAP electrodes for treating vaginal laxity, lubrication issues, and atrophy.



Essence Kit

Includes three reusable CAP intracavitary electrodes with an anti-rotation handle for enhanced pelvic floor treatments.



Single-Use Conductive Intimate Cream

Featuring advanced active ingredients for enhanced comfort and safety.



Longitudinal Electrodes

Designed for deep and superficial pelvic floor therapy with RES and CAP options.

Bibliography

For information on the bibliographic references, please visit the page:

<https://indiba.com/scientific-literature>

If you have any questions, you can contact our medical department at medical@indiba.com

Why Choose INDIBA?

INDIBA's 448 kHz Technology is the gold standard in radiofrequency therapy, trusted by clinics, elite athletes, and rehab professionals worldwide. With 40+ years of research and 150+ clinical studies, it delivers proven results in rehabilitation and intimate health.



Clinically Proven & Expert-Trusted
>80 medical conditions, used globally by top professionals.



Versatile & Effective
Ideal for acute and chronic conditions, ensuring effective results across a wide range of musculoskeletal pathologies.



Safe & Reliable
Engineered for daily clinical use with premium quality.



Education & Support
15+ specialized courses, an online academy, and exclusive post-sales resources via myINDIBA App.



Comfortable & Painless
Ensures a seamless patient experience with measurable results.